


Implant Overdenture System (IOS)

A versatile and innovative solution for long-lasting performance

Implant Overdenture




EACH SET COMES WITH ABUTMENT AND MALE PROCESSING UNIT

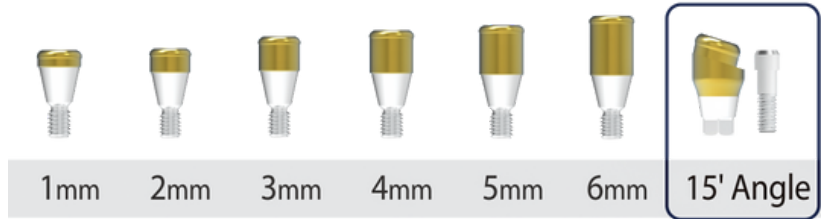

Compatible with most of the major implant systems at tissue, bone, and external levels

- Plastic carrier for the abutment makes it easier, safer to handle and check alignment in the mouse.

- Low vertical height
- Dual Retention - outside and inside retention ensures the longest lasting performance
- Easy Insertion - accurate alignment using self locating design

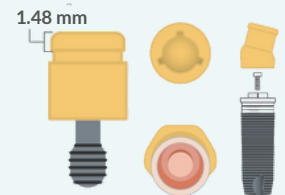


- Dual Retention provides greater resiliency and superior performance








- Resilient attachment with the superior mobility
- 15° Angled abutment can be used with non-parallel implants 1.5mm and 3 mm cuff heights available

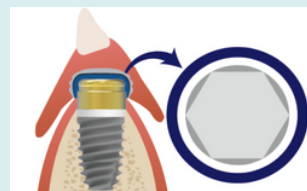

- Advanced metal housing
Minimized hexagonal design prevents rotation and widens surface area of housing for better grip in the denture



Tools and Components

Image	Description	Pkg
	Male Processing Set	Pack of 2 Pack of 10
	Insert Caps - white: 4 lbs retention - Green cap : Angled, 3.2 lbs retention - Pink cap: 2.4 lbs retention - Orange cap: light angled, 1.8 lbs retention - Blue cap: 1.2 lbs retention - Red cap: extra light angled, 0.8 lbs retention - Yellow cap: zero retention Metal housing with inner black plastic cap	Pack of 4 Pack of 10
	Metal Housing	Pack of 4 Pack of 10
	Impression Coping	Pack of 4 Pack of 10
	Lab Analog	Pack of 4 Pack of 10

Optimax Tool (item code: LKW)

Removal End


Turn the body clockwise to spread the removal tip for better grip on the old or worn out caps



Insert End

Combo Torque Wrench Tip (item code: TWT)

Jaw Model

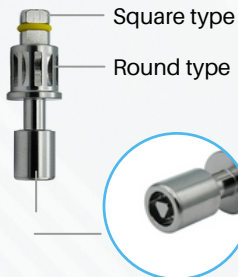

Square type

Round type

Dual head connection (Square & Round) to fit most of the implant system

Round guide makes it easy to attach the driver into the abutment in the mouth